

Ligate & Electroporate 2-4Kb + >4Kb cDNA

Purpose: Make sure the larger sized fractions are okay & determine their sizes.

## Ligation


tube	vector	insert	H <sub>2</sub> O	45°C 5'	19ul 10X ligase buffer 1ul T4 DNA ligase
1	pcDNA1 NotI/BstXI gel purified <span style="border: 1px solid black;">X099</span> 1ul = 24ng	—	8ul	✓	1ul
2	"	cDNA 2-4Kb <span style="border: 1px solid black;">X095</span> 1.5ul = 6ng	6.5ul	✓	1ul
3	"	cDNA >4Kb <span style="border: 1px solid black;">X095</span> 6ul = 6ng	7ul	✓	1ul

15°C overnight

Attachment A

X099 cont'd

## CREATION OF EXPRESSION LIBRARY POOLS

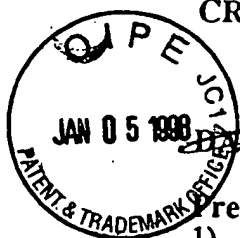

## Prepare bacteria

- 1) ~~Streak out MC1061/P3 onto LB/Kanamycin (15 µg/ml) plate~~

## Prepare reagents and everything else

- 2) Make LB/Amp (15 µg/ml)/Tet (8 µg/ml) plates (15 cm). Make 1 liter total.
- 3) Sterilize:

eppendorfs (1/electroporation)

ddH<sub>2</sub>O (500 ml)

10% glycerol (500 ml)

pasteur pipets (long ones)

1 M MgCl<sub>2</sub>1 M MgSO<sub>4</sub>

1 M glucose

- 4) Put in cold room:
  - all the above sterilized stuff except medium stuff
  - electroporation cuvettes (0.2 cm gap)
  - pipet tips (yellows/blues?)
  - also sign up centrifuges

~~DAY 2:~~

## Bacteria

- 8) Start overnight cultures of MC1061/P3 from A. Pearsons plate of MC1061/P3 from Lodiish lab  
Pick ~~at least 2 colonies~~ into 3 mls each of LB + Kanamycin (15 µg/ml)  
Also streak each onto half a plate of LB/Amp/Tet (only undesirable revertents should grow)

~~DAY 3:~~

## Get bacteria ready for electroporation

- 9) Put 3 ml of starter culture into 250 ml LB + Kanamycin (15 µg/ml)  
Grow until 0.5 - 0.7 O.D.

time	O.D. <sub>600</sub>
7:25	-
10:00	0.320
10:20	0.497
10:35	0.600 ⇒ on ice

## Clean up ligation

- 10) Add TE pH 8.0 to 50 µl
- 11) Add 50 µl Phenol/Chloroform/Isoamly alcohol  
Vortex, Spin and recover top aqueous layer
- 12) Add 50 µl TE pH 8.0 to organic layer to backextract  
Vortex, Spin and recover top aqueous layer and add to previous aq. layer (total = 100 µl)
- 13) Add
  - 50 µl 1X LPA ✓
  - 10 µl 3M NaOAc ✓
  - 20-25 µl 100% ethanol
- 14) Put at -80°C 30 min 8:30
- 15) Spin down at 4°C, remove supe and air-dry (don't dry completely) *rinse w/ 70% EtOH*
- 16) Just before ready to use, resuspend in 2 µl TE (sterile) *Keep on ice. tube dried completely*

Use 2 µl/electroporation  
Freeze rest in cDNA box -20°C.

Get bacteria ready for electroporation (everything on ice!)

- 17) Put culture into ice water to chill 15 min (swirl occasionally) 10:35 - 10:50
- 18) Spin down in 1 disposable conical tube, 4°C, 15 min, 4000 rpm (2600 xg), 10:55 - 11:10
- 19) Decant most but not all liquid (leave equal volume liquid as in pellet). Add 5 ml sterile water and resuspend gently with pipet.
- 20) Add 250 ml ice-cold ddH<sub>2</sub>O (sterile), spin 15 min, 4°C, 4000 rpm
- 21) Repeat steps 19-20 but spin 20 min.
- 22) Pour off as much supe as possible (you'll lose some bugs), add 10% glycerol to 12 mls, gently resuspend cells and spin 8,000 rpm 30 min 4°C in SS-34 (in Falcon 2059 tube)
- 23) Pour off supe getting rid of almost all liquid (you'll lose some cells). You want it thick. Resuspend in 100 µl 10% glycerol (you want it thick)

Used 200 µl, had 100 µl left over  
Next time resuspend in 50 µl? or  
don't add any liquid at all? yes  
should get only 3 electro/250 mls

During spin periods set up for electroporation

- 24) Make SOC from SOB
- 25) Put electroporator chamber on ice
- 26) Connect pulse controller to gene pulser (connect in front the red to red and black to black). The cuvette holder should then be connected to the pulse controller.
- 27) Set to:  
200 ohms  
25 µF  
2.5 kV
- 28) Get everything else ready (Falcon 2059 with 1 ml SOC each, pasteurs, tips, etc)

### Electroporation

- Always do controls: water only (neg. control) and uncut vector (positive control)
- 30) Swirl bacteria with sterile yellow tip. Pipet up 40 µl bacteria to tube #1 on ice. Pipet up and down avoiding generation of bubbles. Let sit 30 sec on ice.
- 31) With fresh tip take up bacteria and put into cuvette as close to bottom as possible without creating bubbles. Quickly shake hard down to bottom (v. important).
- 32) Take off cap, put in electroporator chamber, pulse
- 33) Quickly remove cuvette and add 1 ml SOC. Resuspend with pasteur pipet and transfer to 15 ml round bottom and incubate shaking at 37°C, 60 min. 1 pm - 2 pm
- 34) Repeat steps - for each electroporation.
- 35) Put LB/Amp/Tet plates into hood to dry.
- 36) Plate out 50 - 100 µl/plate to test for electroporation efficiency.  
Use 1:100 of positive control  
1:5 of ligation mix  
undiluted neg. control  
Grow overnight 37°C. Store electroporated bacterial cultures at 4°C up to one week.
- 37) Count colonies.

electro #	DNA	tau	amt plated	colonies	effic. (col/ug)	for 5000 col per plate vol of undi
1	Lig #1 - 2 µl of 8 µl	4.5	10 µl + 90 µl LB	70	$1.2 \times 10^6$	-
2	#2	4.5	↓	2072	$3.5 \times 10^7$	24 µl
3	#3	4.5	↓	1300	$2.2 \times 10^7$	38 µl
4	- (2 µl TE)	4.5	100 µl	13	-	-
5	pcDNA1 + 2.5 Kb insert (1 µl = 0.5 mg) 2 µl	4.5	1 µl + 99 µl LB	553	$5.6 \times 10^8$	-

electro

- #2 pcDNA1 + 2-4 Kb 24 µl/plate ⇒ 42 plates = 210,000 clones
- #3 pcDNA1 + >4 Kb 38 µl/plate ⇒ 26 plates = 130,000 clones

assume 75% loss  
48 mls later  
30 µl/plate ~ 32 plal  
48 µl/plate ~ 20 pla


X099 cont'd

## PCR of colonies to check cDNA sizes

19  $\mu$ l 10X PCR buffer  
 47.5  $\mu$ l T7 primer  
 47.5  $\mu$ l SP6 primer  
 0.95  $\mu$ l Tag polymerase  
 1.52 mix of dNTPs (each) (25mM)  
 73.5  $\mu$ l H<sub>2</sub>O  
 190  $\mu$ l total

3  $\mu$ l dATP  
 3  $\mu$ l dCTP  
 3  $\mu$ l dTTP  
 3  $\mu$ l dGTP  
 12  $\mu$ l dNTPs (freeze in PCR box -20°C)

Aliquot 15  $\mu$ l PCR oil/tube  
 Aliquot 10  $\mu$ l above stock sol'n/tube  
 Flame straight needle, poke colony, then into PCR tube  
 PCR 94°C 30 sec  $\rightarrow$  50°C 30 sec  $\rightarrow$  72°C 2 min  $\times$  35 cycle  
 Add 1  $\mu$ l 10X blue juice  
 Load 6  $\mu$ l/lane onto 0.9% Seakem GT6 agarose minigel

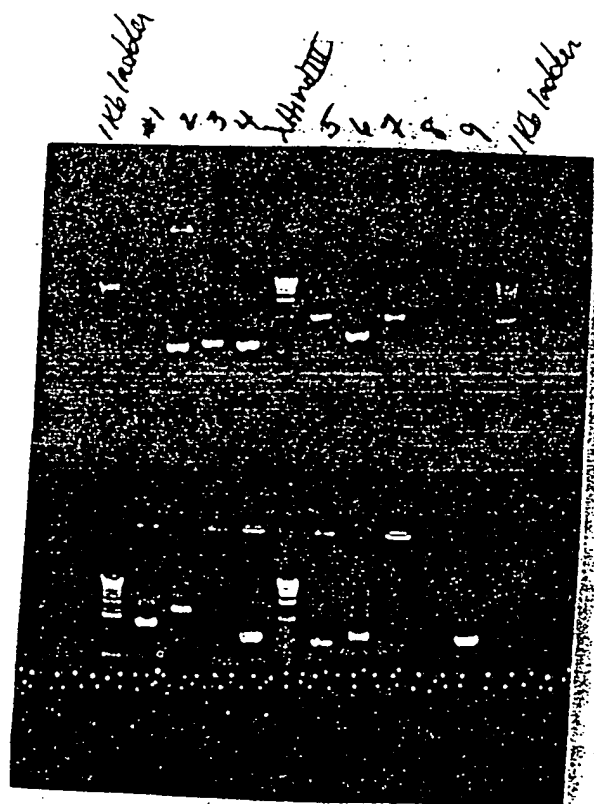

colony	DNA size	minirep
1	-	2.5
2	0.7	
3	0.8	
4	0.7	
5	2.0	
6	1.0	
7	1.9	
8	-	2.1
9	-	2.8
10	1.2	
11	2.4	
12	-	0.8
13	0.9	
14	0.8	
15	1.0	
16	-	3.0
17	-	-
18	0.9	

avg size = 1.6  
 med size = 1.9

avg. size = 1.4  
 med size = 0.95

X099 cont'd

## Mini-preps of clones that did not PCR

Method: Maniatis

- Changes: ① spun twice to get rid of white particulate matter after adding soln III  
 ② Phenol/sevage extracted  
 ③ Resuspended in 25  $\mu$ l TE pH 8

original  
cDNA  
fraction Colonies


2-4 Kb {  
8  
9

> 4 Kb {  
12  
16  
17

Digest	NEB2	10x	10 ng/ml		
plasmid	10x buffer	AcBSA	<del>200 ng/<math>\mu</math>l</del>	H <sub>2</sub> O	No+I+HindIII
<del>2 <math>\mu</math>l</del>	7 $\mu$ l	7 $\mu$ l	0.1 $\mu$ l	25 $\mu$ l	50:50 mix
	+ 7 $\mu$ l				

2  $\mu$ l each + 8  $\mu$ l above rxn mix

37°C 2 hrs (3:45-5:45)

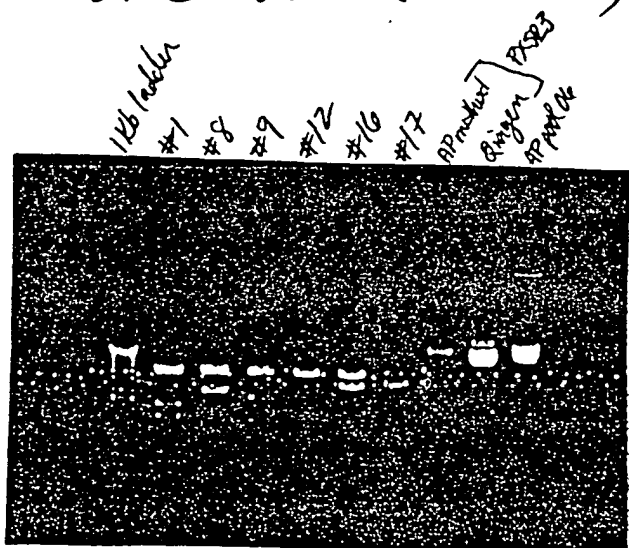

colony	size
#1	1.1 + 1.4 = 2.5
8	2.1
9	1.0 + 1.8 = 2.8
12	0.8
16	3.0
17	-


PX5R23 (1  $\mu$ l each) -  
 Not sure if state of DNA is  
 very different or loading  
 made a difference in the  
 way it ran

## Plasmid midiprep for cDNA library

preps: B46 - B53


## Day 1

1. Scrape 150 mm plate with 5 mls LB. Transfer to Falcon 2059 15 ml tube on ice.
  2. Add another 3 mls LB to plate and scrape again.
  3. Take 400  $\mu$ l, put into freezer vial, add 100  $\mu$ l glycerol and freeze at -150C.
  4. Spin rest in SS-34, 9000 rpm, 2 min 4C.
  5. Dry pellet as much as possible.
  6. Resuspend pellet in 500  $\mu$ l ice cold solution I by vigorous vortexing.
  7. Add 1 ml fresh solution II (0.2 N NaOH, 1% SDS)  
for 100 ml: 1 ml 2 N NaOH  
0.5 ml 20% SDS  
8.5 mls ddH<sub>2</sub>O  
Swirl gently until clear. Do not vortex. Leave on ice 10 min.
  8. Add 750  $\mu$ l solution III (ice-cold). Close tube and mix contents by shaking vigorously several times. Store on ice 5'. A flocculent white precipitate should form.
  9. Centrifuge 15', 4°C, 9000 rpm.
  10. Recover supe and add equal volume of phenol:chloroform. Mix by vortexing.
  11. Spin 9000 rpm, 5'.
  12. Add 2 volumes, ethanol r.t., vortex, let stand 5'.  
Spin in SS-34 rotor for 15' 9,000 rpm.
  13. Carefully aspirate off all liquid (want it dry) and rinse once in cold 70% EtOH. Aspirate again getting any extraneous drops. Airdry 10' minutes.
  14. Redissolve in 100  $\mu$ l TE plus DNase-free RNase (20  $\mu$ g/ml). Vortex briefly. ~~Incubate 37°C, 30'.~~  
Transfer to sterile eppendorf.
  15. Quantitate by dilutions onto EtBr plate.
- Store at 4°C O.N.  
Incubate 37°C 2 1/2 hrs  
Freeze

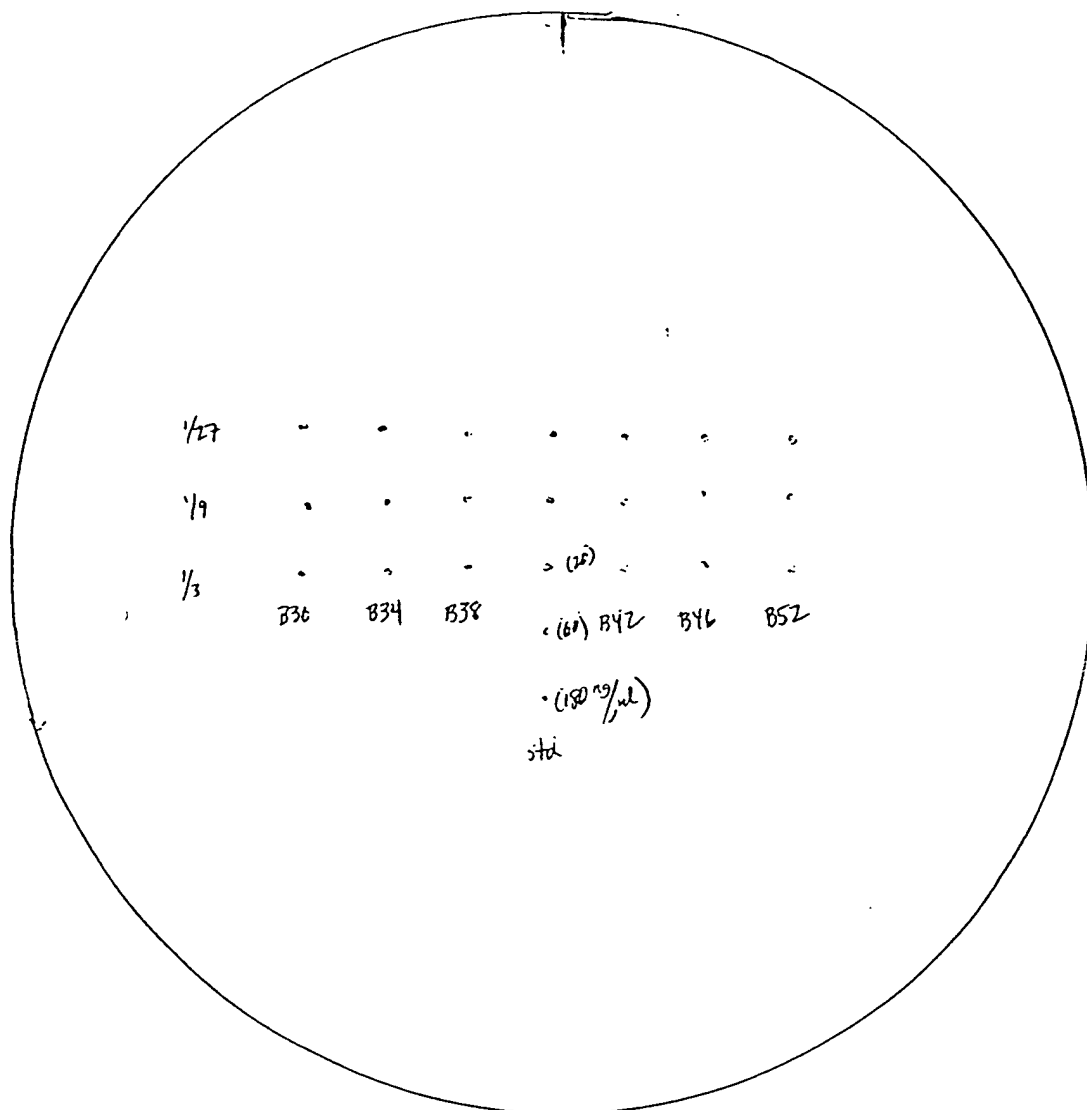

Note:

Photograph is  
mirror image of  
schematic above

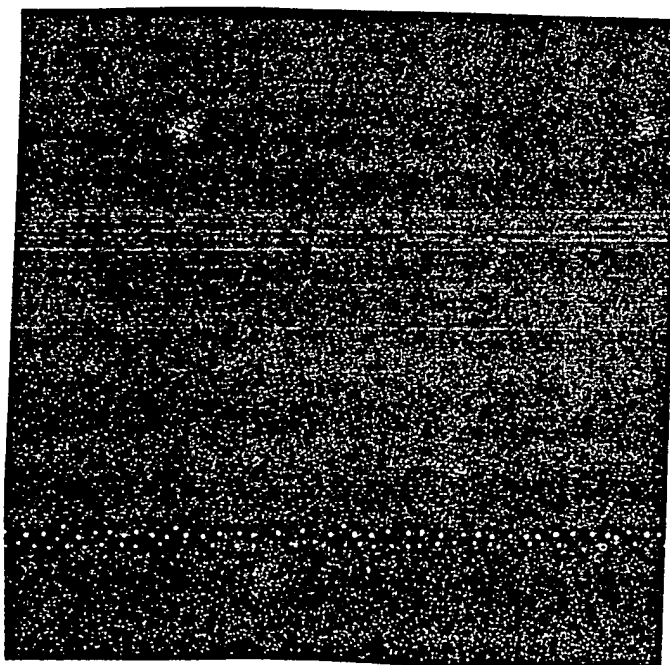

Slight variation  
in preps. I think  
5  $\mu$ l/transfection  
will be good for  
all.

X112

Screen cDNA pools B63-B76, redo B47

## DEAE dextran transfections of COS M6 cells

## materials:

1. 35 mm dishes.
2. DMEM with 10% FBS.
3. Chloroquine (40 mM in CMF PBS, sterile filtered)
4. DNA
5. CMF PBS
6. DEAE-dextran (10 mg/ml in CMF PBS (autoclaved))
7. DMSO
8. cPBS
9. sterile tips

## method:

## day 0 (set up cells)

Set COS M6 cells in 35 mm dishes at 300,000 cells/dish in 2 ml DMEM with 10% FBS

## day 1 (transfect)

1. In sterile eppendorfs prepare for each dish add (in order):

a) DNA - 500 ng/dish

b) add CMF PBS to 190  $\mu$ l, vortexc) 10  $\mu$ l of 10 mg/ml DEAE-dextran, vortex well (7 seconds)

tube #	plates	DNA	CMF PBS	10 mg/ml DEAE-dextran	positives	maybe
1		cDNA pool B63	10 $\mu$ l	180 $\mu$ l	1	2
2		B64			4	7
3		B65			1	0
4		B66			1	0
5		B67			2	0
6		B68			0	1
7		B69			0	4
8		B70			2	3
9		B71			2	4
10		B72			0	1
11		B73			4	0
12		B74			12	4
13		B75			2	2
14		B76			2	5
15		B47 - redo from X111			9	4
16		B47 - redo from X111			12	2
17		pCDNA1	0.39		0	2
18		1:5000	9.2 $\mu$ l		396	8

2. Rinse cells with 2 mls CMF PBS (37C). Aspirate PBS.
3. Add transfection cocktail from step 1 and distribute evenly. Incubate 30 minutes, distributing liquid every 10 min. 9:50 - 10:20
4. Add 2 ml DMEM 10% FBS + 80  $\mu$ M chloroquine and incubate 37C 2.5 hrs. 10:20 - 12:50
5. Aspirate off medium and replace with 1 ml 10% DMSO in DMEM 10% FBS for 2.5 min.
6. Aspirate off and wash once with 2 mls cPBS.
7. Refeed with 2 ml warm DMEM 10% FBS/dish. Incubate overnight.

## Day 2

Refeed cells with 2 ml DMEM 10% FBS + 1 mM Nabutyrate pH 7.3.

## Day 3

Refeed cells with 0.75 ml DMEM 10% FBS + 1 mM Nabutyrate pH 7.3 + 3  $\mu$ g/ml DiI-AcLDL for 5 hrs.14.25 ml DMEM 10% + 14.25  $\mu$ l 1M Nabuty + 158  $\mu$ l DiI-AcLDL prep #48 (0.27 mg/ml)



plate	DNA	positives	maybes
(1)	B47 pool	7 (10)	5 (6)
2	B48	3	4
3	B49	1	0
4	B50	2	6
5	B51	0	1
6	B52	0	0
7	B53	2	2
8	B54	2	5
9	B55	0	0
10	B56	2	4
11	B57	2	5
12	B58	3	8
13	B59	1	4
14	B60	3	2
(15)	B61	6 (4)	2 (3)
16	B62	1	2
17	pcDNA1	1	2
18	1:5000	430	-

Positives are scored if cells are punctate  
parentheses are recounts

X120

Create subpools of 15 colonies of B47.1.8

Purpose: Reduce pool size to approx 15 colonies to narrow the search for the MAC26-1 receptor.

Transform competent MC1061/P3 (Q.G. purple dot) usual procedure but didn't incubate on ice 30'; just heat-shocked 37° 5' right away. Still worked

Plated 5  $\mu$ l - count 930 (B47.1.8)  
4 (no DNA)

Took 3.2  $\mu$ l transformed bugs + 1.9 ml LB  
plated 50  $\mu$ l/plate

Counted plates  
colonies

B47.1.8.1 -	19
.2 -	15
.3 -	19
.4 -	23
.5 -	15
.6 -	13
.7 -	17
.8 -	21
.9 -	15
.10 -	33
.11 -	18
.12 -	17
.13 -	22
.14 -	26
.15 -	11
.16 -	16
.17 -	17
.18 -	15

.19 -	25
.20 -	14
.21 -	19
.22 -	16
.23 -	23
.24 -	23
.25 -	25
.26 -	16
.27 -	19
.28 -	23
.29 -	26
.30 -	22
.31 -	19
.32 -	18
.33 -	17
.34 -	20
.35 -	19
.36 -	31

1st 24 pools  
to 452  
plate 18.8

all 36 pools  
717  
19.9

X120 cont'd

Plasmid midiprep for cDNA library

preps: B47.1.8.1 - B47.1.8.24

Day 1

1. Scrape 100 mm plate with 2 mls LB. Transfer to Falcon 2059 15 ml tube on ice.
2. Add another 2 mls LB to plate and scrape again.
3. Spin in SS-34, 9000 rpm, 2 min 4C.
5. Dry pellet as much as possible.
6. Resuspend pellet in 300  $\mu$ l ice cold solution I by vigorous vortexing.
7. Add 0.6 ml fresh solution II (0.2 N NaOH, 1% SDS)  
for 100 ml: 1 ml 2 N NaOH  
0.5 ml 20% SDS  
8.5 mls ddH<sub>2</sub>O
8. Swirl gently until clear. Do not vortex. Leave on ice 10 min.
8. Add 450  $\mu$ l solution III (ice-cold). Close tube and mix contents by shaking vigorously several times. Store on ice 5'. A flocculent white precipitate should form.
9. Centrifuge 15', 4°C, 9000 rpm.
10. Recover supe and add equal volume (1.2 ml) of phenol:chloroform. Mix by vortexing.
11. Spin 9000 rpm, 5'.
12. Add 2 volumes (2.5) ethanol r.t., vortex, let stand 5'. overnight (1-16) or 6 hrs (17-24)  
Spin in SS-34 rotor for 15' 9,000 rpm.
13. Carefully aspirate off all liquid (want it dry) and rinse once in cold 70% EtOH. Aspirate again getting any extraneous drops. Airdry 10' minutes.
14. Redissolve in 50  $\mu$ l TE plus DNase-free RNase (20  $\mu$ g/ml). Vortex briefly. Incubate 37C, 2 hr. Transfer to sterile eppendorf.
15. Quantitate by dilutions onto EtBr plate.

XIZI

Screen subpools B47.1.8.1 - B47.1.8.24;  
compare CD36 with B47.1.8

DEAE dextran transfections of COS M6 cells

## method:

day 0 (set up cells)

Set COS M6 cells in 35 mm dishes at 300,000 cells/dish in 2 ml DMEM with 10% FBS

day 1 (transfect)

1. In sterile eppendorfs prepare for each dish add (in order):

a) DNA - 500 ng/dish

b) add CMF PBS to 190  $\mu$ l, vortexc) 10  $\mu$ l of 10 mg/ml DEAE-dextran, vortex well (7 seconds)

tube	# plates		CMF PBS	10 mg/ml DEAE-dextran	Results	
1	847.1.8.1		10 $\mu$ l	180 $\mu$ l	10 $\mu$ l	-
2	2					-
3	3					-
4	4					+
5	5					?
6	6					-
7	7					-
8	8					-
9	9					-
10	10					-
11	11					-
12	12					-
13	13					-
14	14					-
15	15					-
16	16					-
17	17					-
18	18					-
19	19					+
20	20					-
21	21					-
22	22					-
23	23					-
24	24					+
25	bSR II (PKSR3)	-	1.5 $\mu$ l	570 $\mu$ l	30 $\mu$ l	v. bright
26		- poly I	X	X	X	nothing
27		mBSA	X	X	X	v. bright
28	CD36	-	1.5 $\mu$ l	525 $\mu$ l	30 $\mu$ l	v. bright
29		- poly I	X	X	X	v. bright
30		mBSA	X	X	X	nothing
31	847.1.8	-	30 $\mu$ l	540 $\mu$ l	30 $\mu$ l	pos. cells
32		- poly I	X	X	X	pos. cells
33		mBSA	X	X	X	nothing?
34	pcDNA1		0.3 $\mu$ l	190 $\mu$ l	10 $\mu$ l	

Pool #4 was brightest &amp; w/ more positive cells.

CD36 binds acetylated LDL & is not inhibited by poly I  
but is inhibited by small amounts of mBSA  
 $\therefore$  CD36 has same properties as MAC 26-1 receptor!

X121 cont'd

2. Rinse cells with 2 mls CMF PBS (37C). Aspirate PBS.
3. Add transfection cocktail from step 1 and distribute evenly. Incubate 30 minutes, distributing liquid every 10 min.
4. Add 2 ml DMEM 10% FBS + 80  $\mu$ M chloroquine and incubate 37C 2.5 hrs. 11:50 - 1:20 pm
5. Aspirate off medium and replace with 1 ml 10% DMSO in DMEM 10% FBS for 2.5 min.
6. Aspirate off and wash once with 2 mls cPBS.
7. Refeed with 2 ml warm DMEM 10% FBS/dish. Incubate overnight.

Day 2

Refeed cells with 2 ml DMEM 10% FBS + 1 mM Nabutyrate pH 7.3.

Day 3

Refeed cells with 0.75 ml DMEM 10% FBS + 1 mM Nabutyrate pH 7.3 + 3  $\mu$ g/ml DiI-AcLDL for 5 hrs.

2.7 mls med + 27  $\mu$ l 1M Nabut + 300  $\mu$ l DiI-AcLDL (#48 0.27  $\mu$ g/ml)

2.25 mls + 225  $\mu$ l poly I (4  $\mu$ g/ml) = 400  $\mu$ g/ml

2.25 mls + 1.35 m-BSA (3.34  $\mu$ g/ml) = 2  $\mu$ g/ml

9:20 - 2:20 pm

# pos cells

4 > 24 > 19

brighter

4 > 24 > 19

X122

Create subpools of 1 colony of B47.1.8 and Screen

Purpose:

Transform competent (Call<sub>2</sub>) MC1061/P3 (Q.6uo purple db)

1. Thaw aliquots (2) of bugs on ice
2. Add 2  $\mu$ l of DNA (or TE) to aliquot  
DNA = pool B47.1.8.4  
TE = neg. control
3. Carefully pipet up & down twice to mix
4. Heat shock 37°C 5 min.
5. Add 200  $\mu$ l LB medium Shake 1 hr 37°C
6. Plate 5  $\mu$ l on 150 mm LB Amp/tet plate

Results: Transf<sup>ormation</sup> worked well. Circled 49 apparently single colonies to be picked. Picked each (done by me Ana Maria Vranceanu) into 3 ml LB A/T. Grew overnight.

Also had Ana Maria Vranceanu repick off same plate & streak onto new plates. (4 indiv/150mm plate)

**X122** cont'd

Plasmid miniprep for cDNA library

preps:

Matrix. 7x7 rows A-F, columns 1-7

14 mini-preps

Day 1

1. Take 200  $\mu$ l culture from each tube in a row or column of 7 and put into eppendorf. Store remainder at 4°C.
2. Spin at 12,000 x g for 30 sec in microfuge.
3. Remove medium by aspiration, leaving bacterial pellet as dry as possible.
4. Resuspend pellet in 100  $\mu$ l ice-cold solution I by vigorous vortexing.
5. Add 200  $\mu$ l fresh solution II (0.2 N NaOH, 1% SDS)  
for 2 ml: 0.2 ml 2 N NaOH  
0.1 ml 20% SDS  
1.7 mls ddH<sub>2</sub>O
6. Swirl gently until clear. Do not vortex. Leave on ice 10 min.  
Add 150  $\mu$ l solution III (ice-cold). Close tube and vortex gently inverted for 5 sec. Store on ice 5'. A flocculent white precipitate should form.
7. Centrifuge 5', 4°C, max speed in microfuge.
8. Recover supe and add equal volume of phenol:chloroform. Mix by vortexing. <sup>400  $\mu$ l</sup>
9. Spin 2' in microfuge. <sup>800  $\mu$ l</sup>
10. Add 2 volumes, ethanol r.t., vortex, let stand 2' at r.t..  
Spin 5', 4°C max speed in microfuge.
11. Carefully aspirate off all liquid (want it dry) and rinse once in cold 70% EtOH. Aspirate again getting any extraneous drops. Airdry 10' minutes.
12. Redissolve in 10  $\mu$ l TE plus DNase-free RNase (20  $\mu$ g/ml). Vortex briefly. Incubate 37°C, 0.5 hr.  
~~Transfer to sterile eppendorf.~~

X122

cont'd

## Screen matrix

## DEAE dextran transfections of COS M6 cells

## materials:

1. 35 mm dishes
2. DMEM with 10% FBS
3. Chloroquine (40 mM in CMF PBS, sterile filtered)
4. DNA
5. CMF PBS
6. DEAE-dextran (10 mg/ml in CMF PBS (autoclaved))
7. DMSO
8. cPBS
9. sterile tips

## method:

day 0 (set up cells) <sup>6-well</sup>  
Set COS M6 cells in 35 mm dishes at 300,000 cells/dish in 2 ml DMEM with 10% FBS <sup>well</sup>

## day 1 (transfect)

1. In sterile eppendorfs prepare for each dish add (in order):
  - a) DNA - 500 ng/dish.
  - b) add CMF PBS to 190  $\mu$ l, vortex.
  - c) 10  $\mu$ l of 10 mg/ml DEAE-dextran, vortex well (7 seconds)

tube #	plates	DNA		CMF PBS	10 mg/ml DEAE-dextran	Results	
1	row A	plasmid prep from matrix		5ul	185ul	10ul	+ 3rd brightest
2	B						-
3	C						-
4	D						+ v. faint few
5	E						-
6	F						+ 2nd brightest
7	G						+ brightest row
8	column 1						+
9	2						+
10	3						-
11	4						+
12	5						+ weak
13	6						-
14	7						+ prob.
15	pcDNA1						+ prob. brighter
16	1:5000 (b5124)		0.3	190ul			-
17			9.3ul	180ul			
18							

2. Rinse cells with 2 mls CMF PBS (37C). Aspirate PBS.
3. Add transfection cocktail from step 1 and distribute evenly. Incubate 30 minutes, distributing liquid every 10 min.
4. Add 2 ml DMEM 10% FBS + 80  $\mu$ M chloroquine and incubate 37C 2.5 hrs.
5. Aspirate off medium and replace with 1 ml 10% DMSO in DMEM 10% FBS for 2.5 min.
6. Aspirate off and wash once with 2 mls cPBS.
7. Refeed with 2 ml warm DMEM 10% FBS/dish. Incubate overnight.

## Day 2

Refeed cells with 2 ml DMEM 10% FBS + 1 mM Nabutyrate pH 7.3.

## Day 3

Refeed cells with 0.75 ml DMEM 10% FBS + 1 mM Nabutyrate pH 7.3 + 3  $\mu$ g/ml DiI-AcLDL for 5 hrs.

13.5 ml

+ 13.5  $\mu$ l 1M Nabut + 150  $\mu$ l DiI-AcLDL #48(0.27  $\mu$ g/ml)

9:20 - 2:20

Results: Lots! of positives. Row G seemed brightest so



# X122 (cont'd) Matrix

	1	2	3	4	5	6	7
A	1	2	3	4	5	6	7
B	8	9	10	11	12	13	14
C	15	16	17	18	19	20	21
D	22	23	24	25	26	27	28
E	29	30	31	32	33	34	35
F	36	37	38	39	40	41	42
G	✓ 43	44	45	46	47	48	✓ 49

As stated before,  
Plated all 49 out on LB Amp/Tet plates  
(before results of this exp<sup>t</sup> were known)

Picked a single colony (hopefully) from  
the plates of # 7, 28, 42, 43, 45, 46, 48, 49

# 42 did not grow, but rest did

Decided to keep things simple so  
mini-preped only # 43, 45, 48, 49 (see X124)

X124 control

## DEAE dextran transfections of COS M6 cells

## materials:

1. 35 mm dishes
2. DMEM with 10% FBS
3. Chloroquine (40 mM in CMF PBS, sterile filtered)
4. DMSO

5. CMF PBS
6. DEAE-dextran (10 mg/ml in CMF PBS (autoclaved))
7. DMSO
8. cPBS
9. sterile tips

## method:

day 0 (set up cells) 6-well well  
Set COS M6 cells in 35 mm dishes at 300,000 cells/dish in 2 ml DMEM with 10% FBS

- day 1 (transfect) - Note: cells were way too heavy (used Coulter counter) Its better to use 1 confluent T75 for 18 wells than use the counter which isn't working.
1. In sterile eppendorfs prepare for each dish add (in order):
    - a) DNA - 500 ng/dish
    - b) add CMF PBS to 190  $\mu$ l, vortex
    - c) 10  $\mu$ l of 10 mg/ml DEAE-dextran, vortex well (7 seconds)

tube #	plates		CMF PBS	10 mg/ml DEAE-dextran	results
1	847.1.8.4.43	- 5 $\mu$ l	185 $\mu$ l	10 $\mu$ l	+
2	.45				-
3	.48				-
4	.49				+ brightest
5	pcDNA1	- 0.5 $\mu$ l	190		-
6	PKSR3	- 0.5 $\mu$ l			+
7					
8					
9					
10					
11					
12					
13					
14					
15					
16					
17					
18					

2. Rinse cells with 2 mls CMF PBS (37C). Aspirate PBS.
3. Add transfection cocktail from step 1 and distribute evenly. Incubate 30 minutes, distributing liquid every 10 min.
4. Add 2 ml DMEM 10% FBS + 80  $\mu$ M chloroquine and incubate 37C 2.5 hrs.
5. Aspirate off medium and replace with 1 ml 10% DMSO in DMEM 10% FBS for 2.5 min.
6. Aspirate off and wash once with 2 mls cPBS.
7. Refeed with 2 ml warm DMEM 10% FBS/dish. Incubate overnight.

Day 2

Refeed cells with 2 ml DMEM 10% FBS + 1 mM Nabutyrate pH 7.3.

Day 3

Refeed cells with 0.75 ml DMEM 10% FBS + 1 mM Nabutyrate pH 7.3 + 3  $\mu$ g/ml DiI-AcLDL for 5 hrs.

6-wells

6.5 wells

0.45  
0.6  
0.5

4.85 ml DMEM

X

0.27 mg/ml

1 ml

1 ml (stock)

B47.1.8.4.49 was chosen as cloned  
MHC2b-1 receptor &  
was renamed phaSR<sup>III</sup>  
(Note: Ana Maria Vincenau  
did several exp'ts to show  
that the plasmid was a  
single one & repeatedly  
gave the expected activity)